
ANIOSYME PRIME

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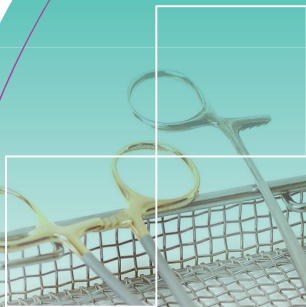UN23670750MCAR-220617
3133002F/17JUN22/S0

ANIOSYME PRIME

EN Cleaning and disinfection of the instrumentation

INDICATIONS: Pre-treatment of surgical instrumentation in foam: during a dry mode of transport or during storage of soiled instruments before machine or manual washing. Maintains an adequate level of humidity to avoid blood and other organic materials drying. **COMPOSITION:**

Quaternary ammonium carbonate, non-ionic surfactants, enzymatic complex (protease, lipase, and amylase), sequestering agent and excipients. **MICROBIOLOGICAL PROPERTIES:** Bactericidal* according to EN 13727/EN 14561 in 15 min. Yeasticidal* according to EN 13624/EN 14562 in 15 min. Tuberculocidal* according to EN 14348/EN 14563 in 30 min. Mycobactericidal* according to EN 14348/EN 14563 in 60 min. Active against enveloped viruses* according to EN 14476/EN 17111 in 15 min. *dirty conditions. **INSTRUCTIONS**

FOR USE: Ready-to-use solution. Beforehand, spray the bottom of the soaking tray with Aniosyme Prime and place instruments collected. Re-cover all the instruments with foam. Cover the tray with the lid. Minimum contact time in the system: 15 minutes. For tuberculocidal activity, allow 30 minutes contact time. Thoroughly rinse treated medical devices with tap water to remove all traces of product and organic soiling before machine or manual cleaning. Not suitable for instruments made of zinc or iron. Scientific dossier (efficacy, compatibility, etc.) available on request.

PRECAUTIONS FOR USE: Contains alcohols, C8-10, ethoxylated. (H318) Causes serious eye damage. (P280) Wear eye/face protective equipment. (P305 + P351 + P338) IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. (P310) Immediately call a POISON CENTER or doctor/physician. Dispose of the product and its container, following all precautions for use. Reserved for professional use. Made in France.

© Keep out of the reach of children. Anti-Poison Center Tel: 210 7793777 (Greece), 1401 (Cyprus).

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