OPASTER'ANIOS

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4P2504922005LEMR-230825
3115360F/25AUG23/50

An Ecolab Company


5L

OPASTER'ANIOS

UFI : C70Y-JR7E
CF0M-NGA8

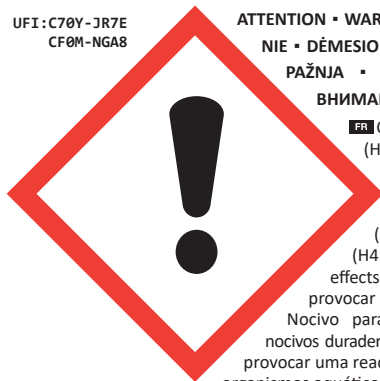ATTENTION • WARNING • ATENCIÓN • ATENÇÃO • OSTRZEŻE-
NIE • DÉMESIO • BRĪDINĀJUMS • UPOZORENJE • POZOR •
PAŽNJA • ΠΡΟΣΟΧΗ • ΔΙΚΚΑΤ • TÄHELEPANU •
ВНИМАНИЕ

FR Contient du 1,2-benzenedicarboxaldehyde.

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EN High level disinfectant for medical devices

INDICATIONS: High level disinfectant for medical devices: surgical and medical equipment, endoscopy equipment and heat-sensitive equipment. **COMPOSITION:** Ortho-Phthalaldehyde solution (0.55%) in a buffered medium, colouring, perfume. **MICROBIOLOGICAL PROPERTIES:** Bactericidal* according to EN 13727/EN 14561 in 5 min. Yeastcidal* according to EN 13624/EN 14562 in 5 min. Fungicidal* according to EN 13624/EN 14562 in 5 min. Tuberculocidal* according to EN 14348/EN 14563 in 5 min. Mycobactericidal* according to EN 14348/EN 14563 in 5 min. Active against enveloped viruses* according to EN 14476/EN 17111 in 5 min. Viricidal* according to EN 14476/EN 17111 in 5 min. Sporicide* according to EN 17126 in 30 min. *=clean conditions. **INSTRUCTIONS FOR USE:** Ready-to-use solution. Do not dilute. Carefully clean and rinse the medical device before soaking it in OPASTER'ANIOS. **Pre-disinfection step:** Clean the medical device with a detergent or pre-disinfectant product. Rinse thoroughly. For endoscopic equipment: brush channels. **Disinfection step:** Pour the solution. Comply with the minimum contact time of 5 minutes. Keep the soaking bath covered during the contact time. For endoscopic equipment: make solution circulate inside channels. Do not add any product into OPASTER'ANIOS. Rinse thoroughly with sterile or filtered water (0.2 µm) to eliminate all traces of the product. For endoscopic equipment: the internal and external parts of the medical device. Dry with a single-use towel. Keep the medical device as aseptically as possible up to use it again. To guarantee the efficacy of the disinfectant, check the Ortho-Phthalaldehyde level with a OPASTER'ANIOS testing strip at the start of the session, when restarting after an interruption and at least once a day. Maximum shelf life of the soaking bath: 14 days except 7 days for France (Instruction no. 220 of 4 July 2016). Scientific dossier (efficacy, compatibility, etc.) available on request. **PRECAUTIONS FOR USE:** (P273) Avoid release to the environment. (P280) Wear protective gloves. (P501) Dispose of unused product and its container as hazardous waste. Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State. Reserved for professional use. Made in France.

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